




MEDISCA[®]
YOUR TRUSTED PARTNER IN COMPOUNDING

PRODUCT NO. 6840

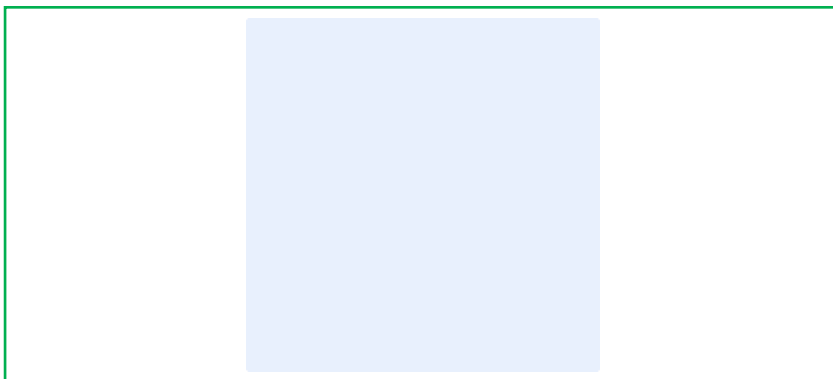

APPLICATOR, TOPI-CLICK, 3 PORT, PINK

35ml, 0.5ml per 2 clicks

DESCRIPTION: The Topi-CLICK 3Port uses the same technology as the Topi-CLICK – but with extra ports for faster and easier CLICK-Dosing of thicker medications. The advantages for the patient are added convenience and ease-of-use when dealing with the thicker and refrigerated medications. The NEW Topi-CLICK 3Port is a smart, effective topical applicator that was designed to overcome the various limitations associated with tubes, tubs, syringes and jars.

FEATURES:

- ✔ Each click moves the plunger approximately 0.25ml of volume space, thereby dispensing 0.25ml of cream on average.
- ✔ The most common number of clicks used for one dose is 2 clicks that move the plunger 0.5ml volume space.
- ✔ Filling with approximately 35ml the Topi-CLICK will provide approximately 64 doses.
- ✔ Around 0.75ml is used for initial priming.
- ✔ Can be labeled like a prescription vial – simplified labeling for lower cost – and easier to read
- ✔ Domed applicator pad for application without dispensing on to hands first.
- ✔ Audible, visual and tactile click registration for user recognition.
- ✔ A “Refill” reminder system helps patients to reorder in a timelier manner.
- ✔ Designed to resist accidental discharge in purse, pocket or luggage.
- ✔ Made of FDA grade materials | BPA and Phthalate free | UV block.
- ✔ One time use only.

TECHNICAL SPECIFICATION:

- ✔ Dispenser Volume: Up to 37 mL
- ✔ Delivery Volume: 5 mL - 30 mL
- ✔ Dispense/CLICK: 0.25 mL
- ✔ Approx. Residual: 2.5 mL
- ✔ Approx. Priming: 0.25 mL or 1 CLICK