



**PRODUCT NO.:6845**



## **Capsule Machine, ProFiller 1100, Fixed Combo w/ 3 sizes: 0,1,3**

**DESCRIPTION:** The fixed combo includes 3 different sizes from 0, 1 and 3, with 1 full system and 2 change parts and tampers to accommodate each size.

Everything you need to fill 3 capsule sizes: Hand-held filler base, locking plate and change parts for 1<sup>st</sup> capsule size and 2 change parts for extra capsule sizes: caps tray, body sheet set (A, B and C), space plate, capsule locker, orienter tray, 50-pin tamper(s), powder spreader, powder tray and powder tray clamps.

### **FEATURES:**

- ✔ Multiple design feature make it easier to use than other 100 hole systems
- ✔ Fill powder using tapping, tamping or vibration (optional)
- ✔ Includes a capsule locker for hand-to-lock powders
- ✔ Ships fully assembled, ready to use in 5 minutes
- ✔ Change capsule sizes in only 2 minutes, change parts available for sizes 000-5
- ✔ Size changes are adjustment and tool free
- ✔ All parts are dishwasher safe and can be dried in a hot air oven or with a portable hair dryer
- ✔ Designed for fast disassembly and visual verification of cleaning
- ✔ High GMP, all metal parts are Inox (316 Stainless Steel)
- ✔ Up to 25% faster than other 100 hole systems based on a patented design and include accessories
- ✔ Inox (Stainless Steel 316) plates for reduced static
- ✔ Tamper size range optimized for higher fill weights and faster filling
- ✔ Unique features include lightweight 2-piece orienter, orienter locator and oversize powder tray and clamps
- ✔ Fully lock capsules after filling, unlike most hand-held fillers – saving you time and frustration


### **APPLICABLE PRODUCTS:**

- ✔ 8314: ProFiller Vibrating Table (100, 3600, 1100)
- ✔ 7352: Capsule Machine, ProFiller 1100, Partial Batch Kit 00-4


### **MATERIAL OF CONSTRUCTION:**

Parts not otherwise labeled are made of Inox (Stainless Steel 316)

**Filler Base with Body Sheet Set** | **Orienter**

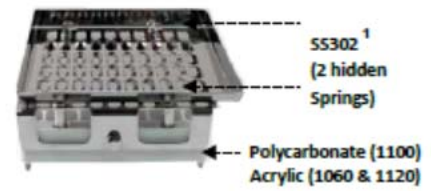

**Caps Tray** | **Locking Plate**



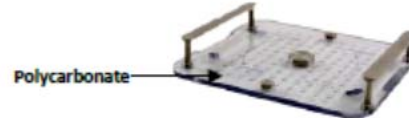


**Accessories**



<sup>1</sup> SS= Stainless Steel (Inox)  
<sup>2</sup> 5-pin Tamper included in item codes P06 HFOAT, P10 HFOAT, and P12 HFOAT

**DIMENSIONS:**

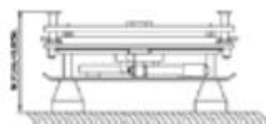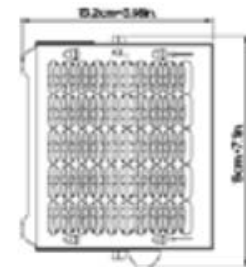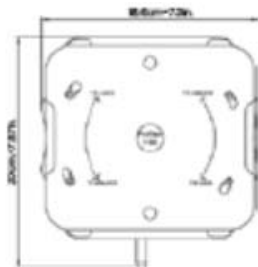

21 x 19 x 10 cm  
8 x 7 x 4 in

18 x 15 x 8 cm  
7 x 6 x 3 in

**Filler**

**Orienter**

USA | CANADA | AUSTRALIA PACIFIC | EUROPE

WWW.MEDISCA.COM

## INSTALLING A CHANGE PART:

### Change Orienter Tray / Orienter Base



1. Loosen 2 thumbscrews and remove Orienter Tray. Place new Orienter Tray and secure with 2 thumbscrews. If changing Orienter Base (changing to/from size 000, 00el, or 5), use new base here.

### Change Body Sheet Set / Spacer Plate



2. Remove Caps Tray. Set Cam Lever to 3 o'clock. Remove 3 Body Sheets.



3. Remove Spacer Plate and replace with new size Spacer Plate. Lubricate Cam and Cam bush on underside of b sheet.

### Change Body Sheet Set / Spacer Plate



4. Place new size Body Sheet a then b, then c on locating pins.



5. Place new size Caps Tray on Filler.



6. Use correct size parts above. You are ready to begin.