



**MEDISCA**<sup>®</sup>  
YOUR TRUSTED PARTNER IN COMPOUNDING

**PRODUCT NO.: 7826**







## **MEDISCA ADAPTER CAP (C, 22mm)**

22-415, LDPE, Green

**DESCRIPTION:** Match the bottle closure-code with the MEDISCA Adapter Cap of your choice for easy connection to bottle tops for simplified oral syringe filling. Please contact a MEDISCA Sales representative for your very own MEDISCA Adaptor-Chek, for quick identification of bottle threads.

### **FEATURES:**

-  Made out of low-density polyethylene (LDPE) with a tethered, re-closable cap
-  Non-DEHP
-  Used on most glass or plastic bottles to draw liquid
-  Latex-free


### **ADDITIONAL INFORMATION:**

#### **Closure Code:**

<b>Letter</b>	<b>Closure Code</b>	<b>mm Code</b>	<b>GPI/SPI Neck Finish</b>
A	18-400	18	400
B	20-400	20	400
<b>C</b>	<b>22-415</b>	<b>22</b>	<b>415</b>
D	24-410	24	410
M	24-415	24	415
E	28-400	28	400
K	28-410	28	410
F	28-415	28	415

## Closure Code Explanation:

- 1) The closure code consists of the cap diameter (mm) and neck finish.
- 2) The neck finishes that are the most common (and that we will work with) are 400, 410, and 415. If you match a bottle neck below with the picture, you will be able to tell the finish.

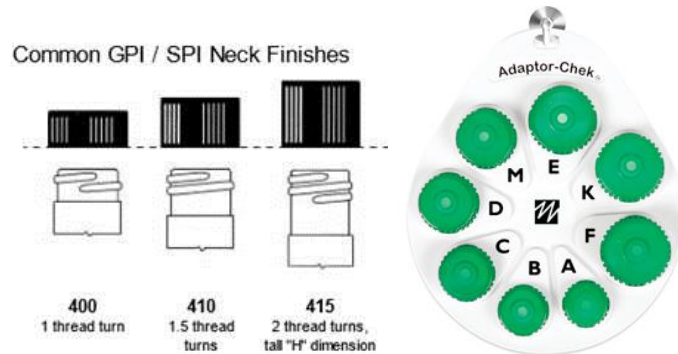

- 3) The closure codes must match, e.g., a 24-410 bottle will only fit a 24-410 cap (adapter-cap D).
- 4) You will see that some Adapter-caps share the same diameter and thus can be easily confused even when using the Adaptor-Chek (seen above). Problems will arise if a D adapter-cap is being used instead of an M, or the like.