



PRODUCT NO.: 7998



PROFILLER 3700 MULTIPLATE MACHINE COMBO A PACKAGE, TORPAC

3 Choice Sizes, Size 00-4 (00, 0, 1, 2, 3 or 4)

DESCRIPTION: Choose 3 different sizes from 00, 0, 1, 2, 3 or 4 & receive 1 full system + 2 change parts.

Everything you need to fill three capsule sizes: Hand-Held Filler Base, Orienter Base, Locking Plate and Change Parts for 1st capsule size / 2 Change Parts for extra capsule sizes: Caps Tray, Body Sheet Set (A, B, C), Spacer Plate, Capsule Locker, Orienter Tray / 300-pin tamper for 1st capsule size (Tampers for additional capsule sizes are sold separately)* / Powder Spreader / Powder Tray + Powder Tray Clamps / Illustrated Manual / Plastic storage box / Empty and Locked Capsule Bins / Set of empty capsule measuring bowls

* Tamper Size Range: 00, 0-1, 2-3, 4

FEATURES:

- Shipped fully assembled. Ready for use in less than 1 hour.
- Easier to use than other 200 or 300 hole systems.
- Dishwasher safe Inox (SS316) powder contact parts. Filler base in lightweight, corrosion resistant anodized aluminum.
- Change capsule sizes in only 5 minutes. Change Parts for Size 000-4 & 00el-2el.
- Orienter Base accommodates capsule sizes 00-4. Size 4 requires Bottom Plate 7180-01 for usage.
- All metal parts Inox (SS316) for reduced static
- Adjustable capsule body height to make over-filling easier & faster. Unique features include lightweight 2-piece orienter (3700), orienter locator and oversize Powder Tray & Push Plate.
- Designed for fast disassembly & visual verification of cleaning.
- Adjustable capsule body height to make over-filling easier & faster.
- Fill weights, USP/EP method within $\pm 1-3\%$ of average. Typically all 300 capsules/batch are within $\pm 4-7\%$ of target. More info at <http://www.profiller.com/3fw>.
- Works with all capsules including gelatin, vegetarian/HPMC and starch capsules.
- Tamper size range optimized for higher fill weights and faster filling.
- Indicator gauges to assist in capsule orientation, separation and locking.
- 3700 fully locks capsules after filling, unlike most hand-held fillers - saving you time and frustration!
- Video tips at www.profiller.com. Well-illustrated manuals in multiple languages.
- Compliance: CE, cGMP (validation package available upon request)

ADDITIONAL INFORMATION:


Warranty: 1 year

Applicable Products:

- 📦 7540: Partial Batch Kit, ProFiller 3600/3700 (Sizes 00-4)
- 📦 7599: Vibrating Table, ProFiller 3700 (110-120V)

Description of Features:

ProFiller 3700 & 3800 Filler

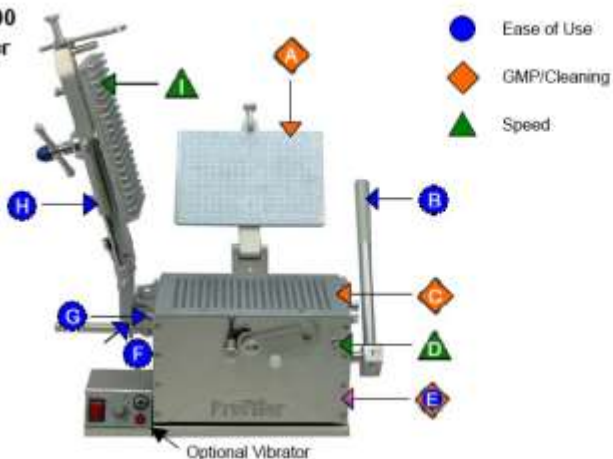


- Ease of Use
- ◆ GMP/Cleaning
- ▲ Speed

◆	SS316 (Inox) Locking Plate. No rubber pad to wear out or clean
●	Rear locking handle for fast, easy movement across front of machine. New internal gear design for easier locking
◆	Dishwasher safe Inox (SS316) powder contact parts including Bottom Plate (not visible)
▲	Adjustable capsule body height for easier & faster over-filling of capsules
●	Base Frame in light-weight anodized aluminum or upgrade to Inox (SS316)
●	Cam position indicator, great for training new operators
●	Over-size thumb screws for easy size change over
●	Dual springs make Tamper easier to lower & raise
▲	Tamper pin size range optimized for faster filling & higher fill weights. Size ranges: 000 / 00 / 0-1 / 2-3 / 4 / 5



◆ Body Sheets fully separate for easy cleaning & visual verification of cleaning

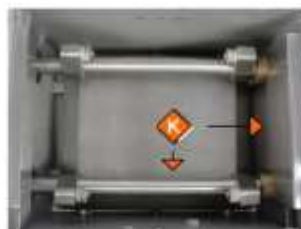


◆ Cover for gear mechanism & Vibrator to minimize powder entry (inside view)



- ▲ Powder Tray is double-depth for faster Spreading
- ▲ Multi-use Push Plate speeds up operations

3700 Orienter

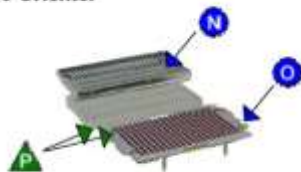

3800 Orienter



●	2-piece Hand-Held Orienter. Shake only the lightweight Orienter Tray	●	Orienter Base common for Size 00-4. Traditional European style Orienter for operators who prefer it
●	Locking indicator confirms capsules are fully locked	●	Guide lines to locate Caps Tray for partial batches
▲	Locating holes in Orienter and tapered locating pins for fast & accurate orientation of capsules into Caps Tray	●	Locking indicator confirms capsules are fully locked

**Adjustment-Free
Size Changeovers**

Shipped Fully Assembled

Dimensions (For Powder Hood Sizing):

