




**PRODUCT NO.: 8027**

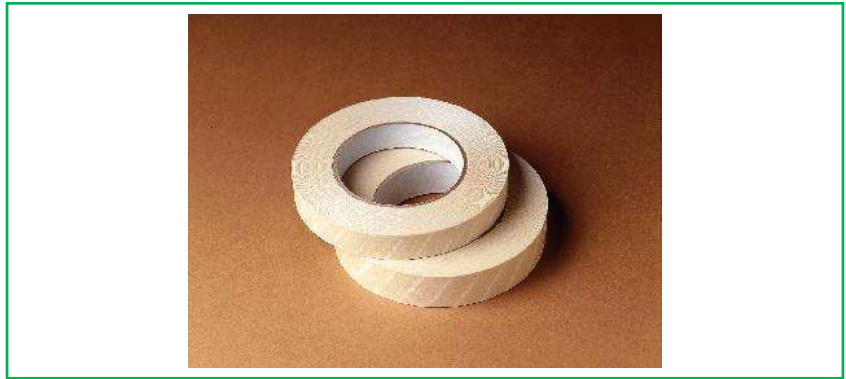

## **AUTOCLAVE INDICATOR TAPE (1" X 60YD)**

2.5cm x 55m

**DESCRIPTION:** Strate-Line's Autoclave Indicator Tape's white stripes turn dark during steam sterilization indicating that the pack has been processed. The strong crepe-embossed paper with latex-free pressure-sensitive adhesive provides secure bonding with fabric and non-woven wraps, paper, metal and glass, yet provides easy removal without traces of residue. High contrast color change ink shows a clear, readily interpreted result, without requiring any guesswork. Marking with permanent ink pen leaves no bleeding on the tape after sterilization.

### **FEATURES:**

- Increased adhesive improves performance during and after sterilization
- Upgraded coating allows for even dispensing and increases ease of handling and application
- Enhanced exterior coating prevents ink transfer and improves chemical indicator appearance for a clearer color change
- Ideal to secure different types of sterilization packs