




MEDISCA® NETWORK INC.  
TECHNICAL SUPPORT SERVICES  
FORMULATION CHEMISTRY DEPARTMENT  
TOLL-FREE: 866-333-7811  
TELEPHONE: 514-905-5096  
FAX: 514-905-5097  
[technicalservices@medisca.net](mailto:technicalservices@medisca.net)

9/23/2020; Page 1

Suggested Formula	Hydroxychloroquine Sulfate 100 mg/5 mL Oral Liquid (Suspension, 100 mL)	FIN	F 003 495v4
-------------------	---	-----	-------------

### SUGGESTED FORMULATION

Ingredient Listing	Qty.	Unit	NDC #	Supplier	Lot Number	Expiry Date
Hydroxychloroquine Sulfate (200 mg) Tablets	10	Units				
Stevia Powder	0.20	g				
Glycerin, USP	5.0	mL				
Tutti Frutti Flavor	0.1	mL				
Chocolate Flavor	0.2	mL				
Medisca Oral Suspend (Suspending Vehicle)	50.0	mL				
Medisca Oral Syrup (Flavored Syrup Vehicle)	q.s. to 100.0	mL				

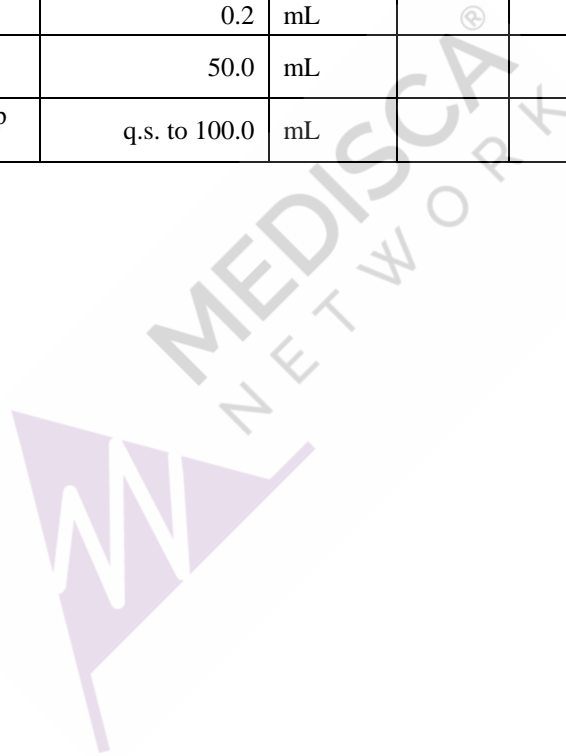



Suggested Formula	Hydroxychloroquine Sulfate 100 mg/5 mL Oral Liquid (Suspension, 100 mL)	FIN	F 003 495v4
-------------------	---	-----	-------------

## SPECIAL PREPARATORY CONSIDERATIONS

### Ingredient-Specific Information

**Light Sensitive** (protect from light whenever possible): *Hydroxychloroquine Sulfate*

**Hygroscopic** (protect from moisture whenever possible): *Stevia Powder, Glycerin*

### Suggested Preparatory Guidelines

Non-Sterile Preparation     Sterile Preparation

Processing Error / Testing Considerations: To account for processing error considerations during preparation, it is suggested to measure an additional **5 to 9%** of the required quantities of ingredients.

Special Instruction: This formula may contain one or more Active Pharmaceutical Ingredients (APIs) that may be classified as hazardous, please refer & verify the current NIOSH list of Antineoplastic and Other Hazardous Drugs in Healthcare Settings. At this time, **General Chapter <800> Hazardous Drugs – Handling in Healthcare Settings** is informational and not compendially applicable unless otherwise specified by regulators and enforcement bodies. For information on the scope, intended applicability, and implementation context for USP General Chapter <800>, see: <https://www.usp.org/compounding/general-chapter-hazardous-drugs-handling-healthcare>.

This formula must be prepared within the appropriate facilities under adequate environmental conditions, following the necessary guidelines and procedures as stated within *USP 795* and *USP 800*, when handling hazardous drugs. Only trained and qualified personnel must prepare this formula.

All required personal protective equipment (hazardous if applicable), such as but not limited to, lab coat, protective sleeves, gloves both inner and outer if applicable, dedicated shoe covers, hairnet, beard cover, eyewear, appropriate face mask, respirator and face shield, etc., where applicable must be worn at all times.

If applicable, follow all required procedures for hazardous drug handling including but not limited to procurement, transport, storage, preparation, dispensing, administration, clean up (spills) & disposal.

If you are a registered 503B facility, please refer to all relevant guidance documents including but not limited to the Code of Federal Regulations (CFR), Guidance for Industry (GFI) and Compliance Policy Guides (CPGs).

This procedure requires the use of very small quantities of ingredients. All calculations and preparation techniques must be verified before dispensing the final product.



Suggested Formula	Hydroxychloroquine Sulfate 100 mg/5 mL Oral Liquid (Suspension, 100 mL)	FIN	F 003 495v4
-------------------	---	-----	-------------

**SUGGESTED PREPARATION (for 100 mL)**

Weigh and / or measure the following ingredients when appropriate:

Ingredient Listing	Qty.	Unit	Multiplication factor (*): _____	Processing Error	Qty. to measure
Hydroxychloroquine Sulfate (200 mg) Tablets §	10	Units			
Stevia Powder §	0.20	g			
Glycerin, USP §	5.0	mL			
Tutti Frutti Flavor	0.1	mL			
Chocolate Flavor	0.2	mL			
Medisca Oral Suspend (Suspending Vehicle)	50.0	mL			
Medisca Oral Syrup (Flavored Syrup Vehicle)	q.s. to 100.0	mL			

§ Weigh / measure just prior to use.

\* Takes into account increased batch size conversions and density conversions, if required.



Suggested Formula	Hydroxychloroquine Sulfate 100 mg/5 mL Oral Liquid (Suspension, 100 mL)	FIN	F 003 495v4
-------------------	---	-----	-------------

**Preparatory Instruction**

**1. Ingredient quantification (determine the actual quantity of Hydroxychloroquine Sulfate (200 mg) tablet powder mix to weigh):**

A. Weigh 11 Hydroxychloroquine Sulfate (200 mg) Tablets. Record the total weight here: \_\_\_\_\_ g

B. Calculate the average weight of powder in each tablet:

Weight of 11 tablets (from Step 1A):	_____ g
DIVIDED BY	
Number of tablets:	11
EQUALS	
Average weight of a single Hydroxychloroquine Sulfate (200 mg) Tablet:	_____ g

C. Calculate the weight of powder equivalent to 10 tablets:

Average weight of a single Hydroxychloroquine Sulfate (200 mg) Tablet (from Step 1B):	_____ g
MULTIPLIED BY	
Number of tablets required:	10
EQUALS	
Weight of powder equivalent to 10 tablets:	_____ g

D. Calculate the weight of powder required *plus* processing error adjustments:

Weight of powder equivalent to 10 tablets (from Step 1C):	_____ g
MULTIPLIED BY	
Processing error adjustments (5% to 9%):	1.05 to 1.09
EQUALS	
<b>Weight of powder required <i>plus</i> processing error adjustments:</b>	_____ g



Suggested Formula	Hydroxychloroquine Sulfate 100 mg/5 mL Oral Liquid (Suspension, 100 mL)	FIN	F 003 495v4
2.	<p><b><u>Powder preparation:</u></b></p> <p>A. Crush and triturate the 11 Hydroxychloroquine Sulfate (200 mg) Tablets into a fine homogeneous powder.</p> <p>B. Sieve and weigh the quantity of Hydroxychloroquine Sulfate (200 mg) tablet powder mix required for the batch (refer to Step 1D) and discard the remaining powder.</p>		
3.	<p><b><u>Powder-liquid preparation:</u></b></p> <p>A. Combine and mix the following ingredients together to form a homogeneous powder blend:</p> <ul style="list-style-type: none"><li>-Hydroxychloroquine Sulfate (200 mg) tablet powder mix (amount weighed in Step 2B)</li><li>-Stevia Powder</li></ul> <p>B. Levigate the homogeneous powder blend (Step 3A) with the Glycerin.</p> <p><u>End result:</u> Homogeneous liquid-like dispersion.</p>		
4.	<p><b><u>Medium Integration:</u></b></p> <p>A. In the given order, sequentially add the following ingredients to the Oral Suspend (Suspending Vehicle):</p> <ul style="list-style-type: none"><li>-Tutti Frutti Flavor</li><li>-Chocolate Flavor</li></ul> <p><u>Specifications:</u> Continuously mix, using high-shear mixing techniques.</p> <p><u>End result:</u> Homogeneous liquid-like dispersion.</p> <p><u>Note:</u> Add the next ingredient, once the previous one has been completely added and dispersed.</p> <p>B. Incrementally add the homogeneous liquid-like dispersion (Step 3B) to the homogeneous liquid-like dispersion (Step 4A).</p> <p><u>Specifications:</u> Continuously mix, using high-shear mixing techniques.</p> <p><u>End result:</u> Homogeneous liquid-like dispersion.</p>		
5.	<p><b><u>Filling to volume:</u></b></p> <p>A. Add Oral Syrup (Flavored Syrup Vehicle) to the mixture (Step 4B) to fill to the required batch size (100.0 mL <i>plus</i> processing error adjustments).</p> <p><u>Specifications:</u> Continuously mix, using high-shear mixing techniques.</p> <p><u>End result:</u> Homogeneous liquid-like dispersion.</p>		



Suggested Formula	Hydroxychloroquine Sulfate 100 mg/5 mL Oral Liquid (Suspension, 100 mL)	FIN	F 003 495v4
-------------------	---	-----	-------------

6.	<p><b><u>Product transfer:</u></b></p> <p>A. Transfer the final product into the specified dispensing container (see “Packaging requirements”).</p> <p><u>Note:</u> Continuously mix the final product during the transfer process into the recommended dispensing containers in order to maintain homogeneity.</p>
----	---

**SUGGESTED PRESENTATION**

Estimated Beyond-Use Date	14 days, refrigerated, as per USP 795.	Packaging Requirements	- Tightly closed, light-resistant dispensing bottle. - To be administered with a metered-dose measuring device.
Auxiliary Labels	1 Use as directed. Do not exceed prescribed dose.	5	Keep refrigerated (2°C – 8°C). Do not freeze.
	2 <b>Shake well before use.</b>	6	Consult your health care practitioner if any prescription or over-the-counter medications are currently being used or are prescribed for future use.
	3 Keep out of reach of children.	7	Cap tightly after use.
	4 May impair mental and/or physical ability. Use care when operating a car or machinery.	8	Protect from light.
Pharmacist Instructions	Add any auxiliary labels specific to the active ingredient to the dispensing container as deemed necessary.		
Patient Instructions	Contact your pharmacist in the event of adverse reactions.		



Suggested Formula	Hydroxychloroquine Sulfate 100 mg/5 mL Oral Liquid (Suspension, 100 mL)	FIN	F 003 495v4
-------------------	---	-----	-------------

## REFERENCES

1.	Glycerin. In: Rowe RC. <i>Handbook of Pharmaceutical Excipients, 4<sup>th</sup> Edition</i> . American Pharmaceutical Association; 2003: 257.
2.	Hydroxychloroquine (Monograph). In: O'Neil MJ. <i>The Merck Index 14<sup>th</sup> Edition</i> . Whitehouse Station, NJ: Merck & Co, Inc.; 2006: Monograph #4818.
3.	Hydroxychloroquine Sulfate (Monograph). <i>United States Pharmacopeia XXXI / National Formulary 26</i> . Rockville, MD. US Pharmacopeial Convention, Inc. 2008.
4.	USP <795>. <i>United States Pharmacopeia XXXI / National Formulary 26</i> . Rockville, MD. US Pharmacopeial Convention, Inc. 2008.
5.	Hydroxychloroquine. Thomson Micromedex. <i>USP DI – Drug Information for the Health Care Professional, 26<sup>th</sup> Edition</i> . Taunton, MA: US Pharmacopeial Convention, Inc; 2006: 1677.

**DISCLAIMER: THIS DOCUMENT IS COPYRIGHT© 2019-2020 MEDISCA PHARMACEUTIQUE INC. MEDISCA NETWORK INC., HEREBY REFERRED TO AS 'THE NETWORK', HAS PROVIDED THE FORMULA AND INSTRUCTIONS ABOVE AS A MODEL FOR EDUCATIONAL PURPOSES ONLY ON THE BASIS OF THE RECOGNIZED COMPENDIA AND TEXTS REFERENCED AT THE END OF THIS DOCUMENT. THE NETWORK TAKES NO RESPONSIBILITY FOR THE VALIDITY, SCHEDULING OR ACCURACY OF THIS INFORMATION OR FOR ITS SAFETY OR EFFECTIVENESS, NOR FOR ANY USE THEREOF, WHICH IS AT THE SOLE RISK OF THE LICENSED PHARMACIST OR OTHER APPROPRIATELY STATE LICENSED PROFESSIONAL. ADJUSTMENTS MAY BE NEEDED TO MEET SPECIFIC PATIENT NEEDS AND IN ACCORDANCE WITH A LICENSED PRESCRIBER'S PRESCRIPTION. THE PHARMACIST OR OTHER APPROPRIATELY STATE LICENSED PROFESSIONAL MUST EMPLOY APPROPRIATE TESTS TO DETERMINE THE STABILITY OF THIS SUGGESTED FORMULA. THE NETWORK CANNOT BE HELD LIABLE TO ANY PERSON OR ENTITY CONCERNING CLAIMS, LOSS, OR DAMAGE CAUSED BY, OR ALLEGED TO BE CAUSED BY, DIRECTLY OR INDIRECTLY, THE USE OR MISUSE OF THE INFORMATION CONTAINED IN THIS SUGGESTED FORMULA. IN ALL CASES IT IS THE RESPONSIBILITY OF THE LICENSED PHARMACIST OR OTHER APPROPRIATELY STATE LICENSED PROFESSIONAL TO KNOW THE LAW, TO COMPOUND ANY FINISHED PRODUCT AND TO DISPENSE THESE PRODUCTS IN ACCORDANCE WITH FEDERAL AND STATE LAW. MEDISCA NETWORK INC. MAKES NO WARRANTIES WITH RESPECT TO INFRINGEMENT OR NON-INFRINGEMENT BY THE FORMULA OF ANY PATENT OR OTHER INTELLECTUAL PROPERTY OF ANY OTHER PARTY, AND IT IS THE RESPONSIBILITY OF THE PHARMACIST OR OTHER APPROPRIATELY STATE LICENSED PROFESSIONAL TO INVESTIGATE AND DETERMINE ANY SUCH ISSUE.**